


Dr Mark's HyGenie® HyGenie Sport



How To Use Your HyGenie Device Find our User Guide Video's at: www.drmarkshygenie.com


1 OPEN

Hold the bottom handle (blue) and device body (clear) in one hand. Turn the top handle (white or red) either direction until it clicks into place. Lift off the top.

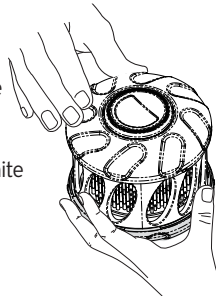


2 PLACEMENT

Nestle your dental appliance into the bottom brush. Paired orthodontic aligners can be placed together, with one facing upwards and one facing downwards, as shown.

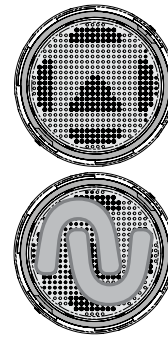


3 CLOSE

Line up the arrows found on the top handle and clear body, as shown. Gently press down until the handle 'sits' into place. Twist it in either direction and it will 'click & lock' closed.

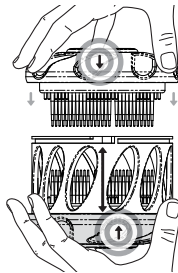


4 DENTAL FRESH

Holding your HyGenie over a basin, squeeze 2-3 generous drops of Dr Mark's DentalFresh® into the device through the perforations and apply a splash of water.

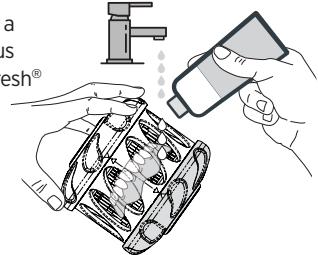


5 WASH

Rotate the top and bottom handles in opposite directions, or, alternate top and bottom washing by holding one handle in place and rotating the opposite handle. Wash for 30-45 seconds or until DentalFresh foam thoroughly covers your appliance.

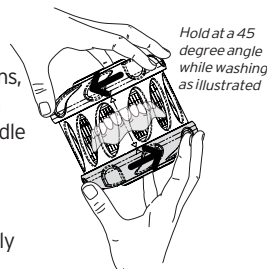


6 RINSE

Rinse the HyGenie under gently running water, rotating the handles in opposite directions, until all DentalFresh foam is rinsed away and water runs clear through the device.

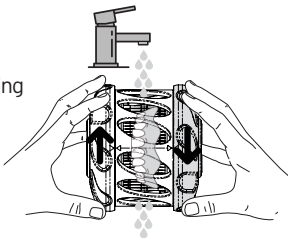

7 WEAR IT, STORE IT, PROTECT IT

Now that your dental appliance is clean, fresh and ready, you can wear it or store it safely in Dr Mark's HyGenie® until you need it next.

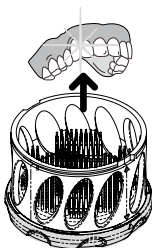

CARE & TIPS Your HyGenie package also doubles as a soaking bath and a handy travel case for 'HyGenie On The Go'. Dr Mark's HyGenie® is self-cleaning and self-drying, just like your toothbrush.

If ever needed, clean your HyGenie brushes by removing both handles, apply DentalFresh onto the brushes, rub together until foamy then rinse. We recommend replacing your Dr Mark's HyGenie® brushes every six months. Order at drmarkshygenie.com

HOW CAN WE HELP?

For more advice on removable dental appliance care and information on Dr Mark's other products, email us at customercare@drmarkshygenie.com

Dr Mark's HyGenie® and HyGenie Sport® are registered trademarks of Dr Mark's HyGenie Pty Ltd.
2.12, 32 Delhi Road, North Ryde, NSW 2113 Australia
ABN 64604680 891

